

Helicobacter Pylori

H. pylori Agtigen Rapid Test Kit Package Insert

A rapid test for the qualitative detection of *Helicobacter pylori* in faeces.
For professional *in vitro* diagnostic use only.

IVD

INTENDED USE

This H. pylori Agtigen Rapid Test Kit is intended for the direct qualitative detection of the presence of H. pylori antigen in patient fecal samples. The test might be used as an aid for detecting patients with acute and chronic gastroenteritis infected with H. pylori.

PRINCIPLE

The H. pylori Agtigen Rapid Test Kit is a sandwich lateral flow chromatographic immunoassay. The test strip consists of: 1) a burgundy colored conjugate pad containing monoclonal anti-H.pylori antibody conjugated with colloid gold (anti-H.P conjugates) and 2) a nitrocellulose membrane strip containing a test band (T band) and a control band (C band). The T band is pre-coated with another monoclonal anti-H.P antibody, and the C band is pre-coated with goat anti-mouse IgG antibody.

When an adequate volume of extracted fecal specimen is dispensed into the sample well of the test cassette, the specimen migrates by capillary action across the cassette. H.P antigen if present in the specimen will bind to the anti-H.P conjugates. The immunocomplex is then captured on the membrane by the pre-coated antibody, forming a burgundy colored T band, indicating a H.P positive test result. Absence of this band suggests that the concentration of H.P in the specimen is below the detectable level, indicating a H.P negative result.

Absence of the T band suggests a negative result. The test contains an internal control (C band) which should exhibit a burgundy colored band of the immunocomplex of goat anti-mouse IgG/mouse IgG-gold conjugate regardless of the color development on the T band. Otherwise, the test result is invalid and the specimen must be retested with another device.

PACKAGING CONTENTS

- 25 H. pylori Agtigen Rapid Test Kit
- 25 Stool collection tubes with buffer
- 25 Plastic droppers for transferring watery stool

STORAGE AND STABILITY

Store the test kit at 2~30°C with a valid period of 24 months. The test kit is stable through the expiration date printed on the sealed pouch.

The test kit must remain in the sealed pouch until use. Use the test kit within 1 hour once the foil pouch is opened.

DO NOT FREEZE or expose the kit over 30°C.

Do not use beyond the expiration date.


SPECIMEN COLLECTION AND PREPARATION

1. Collect a random stool sample in a clean, dry receptacle.

Solid stool samples


2. Open the stool collection device by unscrewing the top and use the collection stick to randomly pierce the stool sample in at least five different sites. Do not scoop stool sample as this may lead to an invalid test result.
3. Ensure stool sample is only in the grooves of the collection stick. Excess stool sample may lead to an invalid test result.
4. Replace the collection stick and tighten securely to close the stool collection device.
5. Shake the stool collection device vigorously.
6. The specimen is now ready for testing, transportation or storage.

Note: Specimens extracted may be stored at 2°C-8°C for up to 3 days. If longer storage is required, freezing at ≤-20°C is recommended.

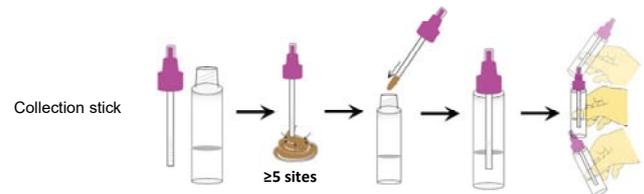

Watery stool samples

5. Collect a random stool sample in a clean, dry receptacle.
6. Open the stool collection device by unscrewing the top.
7. Fill the plastic dropper with the sample; dispense 2 drops (70-85 µL) into the stool collection device.
8. Replace the collection stick and tighten securely to close the stool collection device.
9. Shake the stool collection device vigorously.

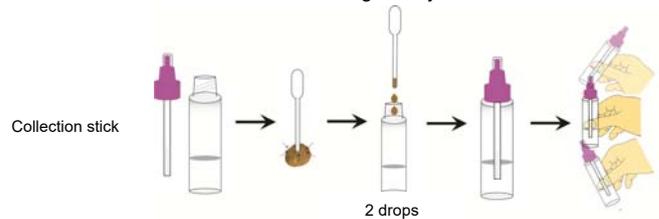


MATERIALS REQUIRED

- Specimen collection container
- Timer
- Positive Control and Negative Control

PROCEDURE

Allow the tests, stool samples and controls to reach to room temperature (15-30°C) prior to testing. Do not open pouches until ready to perform the assay.

1. When ready to test, open the pouch at the notch and remove the test device. Place the test device on a clean, flat surface.
2. Shake the stool collection device vigorously to ensure a homogenous liquid suspension.
3. Position the stool collection device upright and twist off the dispenser cap. Holding the stool collection device vertically, dispense 2 drops of the solution into the sample well of the test device. Do not overload sample.

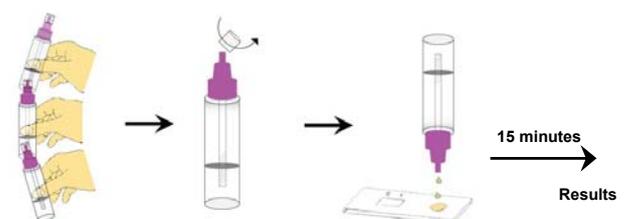

4. Set up the timer.

5. Results can be read in 15 minutes after adding the specimen. Positive results can be visible in as short as 1 minute.

Don't read results after 20 minutes. To avoid confusion, discard the test device after interpreting the result.

INTERPRETATION OF RESULTS

NEGATIVE RESULT: If only the C band is developed, the test indicates that no detectable *H. Pylori antigen* is present in the specimen. The result is negative.

POSITIVE RESULT: If both C and T bands are developed, the test indicates for the presence of *H. Pylori antigen* in the specimen. The result is positive.

Samples with positive results should be interpreted in conjunction with other testing procedure and clinical findings before a diagnostic decision is made.

INVALID: If no C line is developed, the assay is invalid regardless of any color development on the T line as indicated below. Repeat the assay with a new test device. **Excess fecal specimen can lead to invalid test results; if this is the cause, re-sample and re-test (see instructions for collection of specimen).**



QUALITY CONTROL

Internal procedural control is included in the test. A Control line appearing in the control region (C) is the internal procedural control.

It confirms sufficient specimen volume and correct procedural technique. A clear background is an internal negative background control. If the test is working properly, the background in the result area should be clear and not interfere with the ability to read the result.

It is recommended that a positive and a negative control be evaluated to verify proper test performance when:

1. A new operator uses the kit.
2. A new test kit is used.
3. A new shipment of kits is used.
4. The temperature during storage of the kit falls outside of 2-30°C.
5. The temperature of the test area falls outside of 15-30°C.

CLINICAL SIGNIFICANCE

Helicobacter pylori is associated with a variety of gastrointestinal diseases including non-ulcer dyspepsia, duodenal and gastric ulcer and active, chronic gastritis^{1,2}. The prevalence of *H. pylori* infection could exceed 90% in patients with signs and symptoms of gastrointestinal diseases. Recent studies indicate an association of *H. Pylori* infection with stomach cancer³.

H. pylori can be transmitted through the ingestion of food or water that is tainted with fecal matter. Antibiotics in combination with bismuth compounds showed to be effective in treating active *H. Pylori* infection.

H. pylori infection is currently detected by invasive testing methods (ie histology, culture) based on endoscopy and biopsy, or non-invasive testing method, such as urea breath test (UBT), serologic antibody test and stool antigen test. UBT test requires expensive lab equipment and radioactive reagent. Serologic antibody tests do not distinguish between currently active infection and past exposure or infection that has been cured. The stool antigen test detects antigen presence in the feces that indicates an active *H. pylori* infection. It can also be used to monitor the effectiveness of treatment and the recurrence of the infection.

The *H. pylori* Agtigen Rapid Test Kit uses a combination of the polyclonal anti *H. pylori* antibody and collid gold conjugated monoclonal anti-*H. Pylori* antibody to specifically detect *H. pylori* antigen present in the infected patient fecal specimen. The test is user friendly, accurate, and the result is available instantly.

ANALYTICAL PERFORMANCE

Clinical Performance

324 fecal samples collected from subjects with symptomatic gastrointestinal disorders and non-gastrointestinal symptoms were tested with *H. pylori* Agtigen Rapid Test Kit and with the UBT as reference test. A comparison of the results for all subjects is shown in the following table:

| | H. pylori Agtigen Rapid Test Kit | | |
|--------------|----------------------------------|----------|-------|
| | Positive | Negative | Total |
| UBT Positive | 118 | 7 | 125 |
| UBT Negative | 0 | 199 | 199 |
| Total | 118 | 206 | 324 |

Relative Sensitivity: 94.4% , Relative Specificity: 100.0%, Overall Agreement: 97.8%

The detection limit for *H. pylori* Agtigen Rapid Test Kit is 5 ng/ml of *H. pylori* lysate. Fecal specimen extractions containing *H. pylori* lysate equal to or greater than 5 ng/ml routinely test positive. Specimens containing *H. pylori* lysate less than 5 ng/ml may also produce a very faint positive line, especially with an assay time extended beyond 15 minutes.

The following experiments were done to validate the sensitivity of the *H. pylori* Agtigen Rapid Test Kit:

Normal fecal specimen extractions were spiked with *H. pylori* lysate to concentrations of 0, 1.25, 2.5, 5, 10, 20 ng/ml. The specimens were run on the *H. pylori* Agtigen Rapid Test Kit.

Results are shown in the table below.

| H. pylori/lysate ng/ml | 0 | 1.25 | 2.5 | 5 | 10 | 20 |
|------------------------|----|------|-----|----|----|----|
| Number of positive | 0 | 0 | 12 | 20 | 20 | 20 |
| Number of negative | 20 | 20 | 8 | 0 | 0 | 0 |

n=20 Relative Sensitivity at 5 ng/ml = 20/20 x 100% = 100%

PRECAUTIONS

- This package insert must be read completely before performing the test. Failure to follow the insert gives inaccurate test results.
- Do not open the sealed pouch, unless ready to conduct the assay.
- Do not use any kit components beyond their stated expiration dates.
- Do not use the components from any other type of test kit as a substitute for the components in this kit.
- Bring all reagents to room temperature (15°C-30°C) before use.
- **Do not scoop stool sample as this may lead to excess fecal specimen that tends to clot the sample pad and interfere with sample migration.**
- Wear protective clothing and disposable gloves while handling the kit reagents and clinical specimens. Wash hands thoroughly after performing the test.
- Users of this test should follow the US CDC Universal Precautions for biosafety.
- Do not smoke, drink or eat in areas where specimens or kit reagents are being handled.
- Extraction buffer contains 0.1% NaN₃. Avoid contact with skin or eyes. Do not ingest.
- Dispose of all specimens and materials used to perform the test as biohazardous waste.
- The test results should be read within 15 minutes after a specimen is applied to the sample well of the device. Reading results after 20 minutes may give erroneous results.
- Do not perform the test in a room with strong air flow, i.e. electric fan or strong air-conditioning.

LIMITATIONS

1. The Assay Procedure and the Interpretation of Assay Result sections must be followed closely when testing for the presence of *H. pylori* antigen in feces. Failure to follow the procedure, particularly the Specimen Collection procedure, may cause inaccurate results.
2. The assay is limited to the qualitative detection of *H. pylori* antigen in human fecal specimen. The intensity of the test line does not have a linear correlation with the antigen titer in the specimen.
3. A negative result for an individual subject indicates the absence of detectable *H. pylori* antigen. However, a negative test result does not preclude the possibility of infection with *H. pylori*.
4. A negative result can occur if the quantity of the *H. pylori* antigen present in the specimen is below the detection limits of the assay or if the antigens that are detected are not present in the fecal sample collected.
5. If symptoms persist and the result from the *H. pylori* Agtigen Rapid Test Kit is negative or non-reactive, it is recommended to re-sample the patient a few days later or test with alternative test methods.
6. The results obtained with this test should only be interpreted in conjunction with other diagnostic procedures and clinical findings.

REFERENCES

1. Marshall BJ, et al, Pyloric campylobacter infection and gastroduodenal disease. Med.J. Aust. 1985, 142: 439-444
2. Lambert IR, Lin SK, and Aranda-Michel. J, helicobacter Pylori Scan. J. Gastroenterol. 1995, 30 suppl 208: 33-46
3. Vans DJ, Evans DG, et al A sensitive and specific serologic test for detection of campylobacter pylori infection. Gastroenterology. 1989, 96:1004
4. Shimoyama T, KatoC, et al. Applicability of a monoclonal antibody-based stool antigen test to evaluate the results of Helicobacter pylori eradication therapy. 2009, May 62(3): 225
5. Krausse R, Muller G, Doniec M. Evaluation of a rapid new stool antigen test for diagnosis of Helicobacter pylori infection in adult 1008, 46(6): 2062
6. Altman E, Fernandez H. et al Analysis of Helicobacter pylori isolates from Chile: Occurrence of selective type I Lewis b antigen expression in lipopolysaccharide. 2008, 57(pt 5): 585

INDEX OF SYMBOLS

| | | | |
|---|---|---|------------------------------------|
|  | Caution |  | In Vitro Diagnostic Medical Device |
|  | Manufacturer |  | Date of Manufacture |
|  | CE Marking |  | Do Not Re-use |
|  | Keep Dry |  | Keep Away From Sunlight |
|  | Batch Code |  | Do Not Use if Package is Damaged |
|  | Catalogue Number |  | Contains Sufficient for <n> Tests |
|  | Use-By Date |  | Temperature Limit |
|  | Authorized representative in the European Community | | |



Manufacturer

Jinan Chuangxingwell Biotech Co., Ltd.
 Add.: Building 1-2B-402, Liando U Valley,
 Gangyuan Sixth Road, Jinan Zhangjin
 Comprehensive Bonded Zone, Gaoxin
 District, Jinan city, Shandong Province,
 China

Tel: +86-0531-88894020

Email: s@1stiot.com

EC REP

Riomavix S.L.

Add.: Calle de Almansa 55, 1D, Madrid 28039 Spain

E-mail: leis@riomavix.com

Tel.: +34 658 396 230

(SRN: ES-AR-000001202)